


Brightfield Healthcare 9000

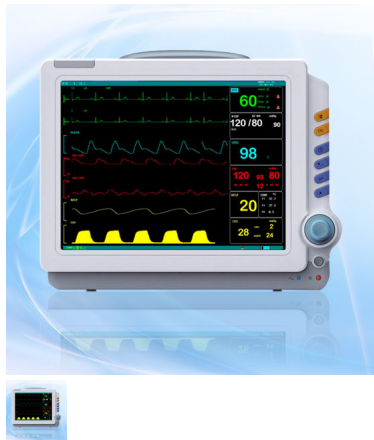

12.1" portable patient monitor

Rating: Not Rated Yet

[Ask a question about this product](#)

Description

12.1" portable patient monitor

Portable patient Monitor is adaptable to adult, pediatric and neonatal usage. It can monitor vital signals as ECG, Respiratory Rate, SpO2, NIBP and TEMP. It integrates parameter measuring modules, display and recorder in one device, featuring in compactness, lightweight and portability. Replaceable built-in battery facilitates transportation of patient. Large high-resolution display provides clear view of 5 waveforms and full monitoring parameters.

Features

- 12.1" color TFT LCD screen, resolution: 800*600
- 5 kinds of interfaces: Standard/Big font/6-lead ECG/Trend graph/OxyCRG
- Maximum 8-channel waveforms display, the waveform color can be defined by user (7 colors)
- 96-hour storage and review of trend graphic gram and table, 400 groups NIBP data, 1800 alarm events
- ST and arrhythmia analysis, pacemaker analysis, drug dose calculation
- Anti-flicker design
- NIBP over pressure protection
- Plug-in rechargeable battery, AC/DC available
- Support wire or wireless network